




1 The polymer core expands when exposed to water, effectively increasing the packing density, and achieving compact embolization, effectively reducing the long-term recurrence rate.

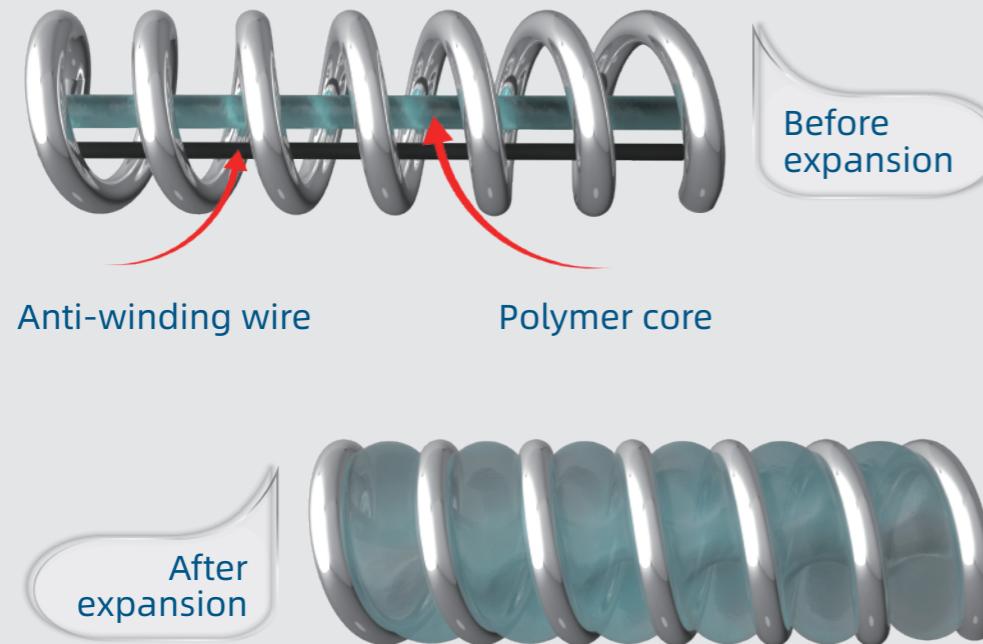

2 The expansion of the polymer core can effectively increase the stability of the coil and avoid the compression of the packed coil under the impact of blood flow, thereby ensuring long-term effective sealing of the aneurysm neck, eliminating dog-ear remnant and effectively reducing the recurrence rate.

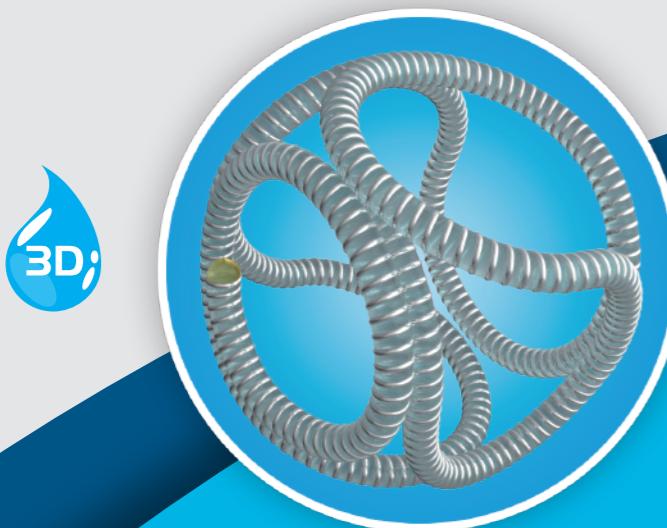

3 Comparing with bare metal coils, Perfiller of the same length has a lower metal occupancy rate, a softer coil, and improved efficacy and safety during use.

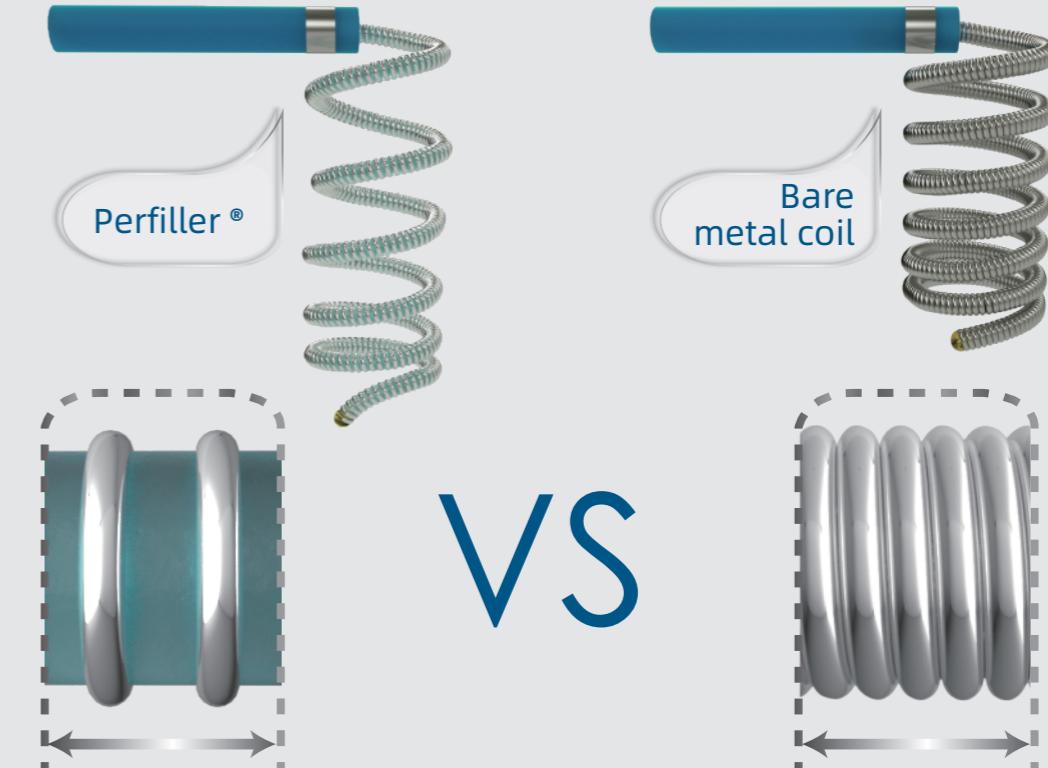

4 Minimum loop diameter is 1mm and a wide range of small size coils provide more options for small aneurysms as finishing coils.



5 The new delivery system has greater pushability and crossability. The torque response and transmission of the pusher are improved. The delivery process is easier and smoother.

This product is indicated for the endovascular embolization of intracranial aneurysms and other neurovascular abnormalities such as arteriovenous malformations and arteriovenous fistulae.

Product Specification

Coil Structure		3D			2D			
Softness	Specification	Loop Dia. (mm)	Length (cm)	Coil Structure	Specification	Loop Dia. (mm)	Length (cm)	Coil Structure
Extra Soft	TJCHYF0102	1	2	Complex	TJCHYL0102	1	2	Helical
	TJCHYF1.502	1.5	2	Complex	TJCHYL1.501	1.5	1	Helical
	TJCHYF1.503	1.5	3	Complex	TJCHYL1.502	1.5	2	Helical
	TJCHYF0202	2	2	Complex	TJCHYL1.504	1.5	4	Helical
	TJCHYF0204	2	4	Complex	TJCHYL0201	2	1	Helical
	TJCHYF0206	2	6	Complex	TJCHYL0202	2	2	Helical
	TJCHYF2.502	2.5	2	Complex	TJCHYL0203	2	3	Helical
	TJCHYF2.504	2.5	4	Complex	TJCHYL0204	2	4	Helical
	TJCHYF2.506	2.5	6	Complex	TJCHYL0206	2	6	Helical
	TJCHYF0304	3	4	Complex	TJCHYL0208	2	8	Helical
	TJCHYF0306	3	6	Complex	TJCHYL2.502	2.5	2	Helical
	TJCHYF0308	3	8	Complex	TJCHYL2.504	2.5	4	Helical
	TJCHYF0310	3	10	Complex	TJCHYL2.506	2.5	6	Helical
	TJCHYF0405	4	5	Complex	TJCHYL0302	3	2	Helical
	TJCHYF0408	4	8	Complex	TJCHYL0303	3	3	Helical
	TJCHYF0410	4	10	Complex	TJCHYL0304	3	4	Helical
Soft	TJCHYF0510	5	10	Complex	TJCHYL0306	3	6	Helical
	TJCHYF0515	5	15	Complex	TJCHYL0308	3	8	Helical
	TJCHYF0612	6	12	Complex	TJCHYL0310	3	10	Helical
	TJCHYF0620	6	20	Complex	TJCHYL0312	3	12	Helical
	TJCHYF0715	7	15	Complex	TJCHYL0404	4	4	Helical
	TJCHYF0720	7	20	Complex	TJCHYL0406	4	6	Helical
	TJCHYF0728	7	28	Complex	TJCHYL0408	4	8	Helical
	TJCHYF0817	8	17	Complex	TJCHYL0410	4	10	Helical
	TJCHYF0830	8	30	Complex	TJCHYL0415	4	15	Helical
	TJCHYF0920	9	20	Complex				
	TJCHYF0930	9	30	Complex				
	TJCHYF1020	10	20	Complex				
	TJCHYF1030	10	30	Complex				
	TJCHYF1120	11	20	Complex				
	TJCHYF1130	11	30	Complex				
	TJCHYF1140	11	40	Complex				
Standard	TJCHYF1225	12	25	Complex				
	TJCHYF1325	13	25	Complex				
	TJCHYF1345	13	45	Complex				
	TJCHYF1430	14	30	Complex				
	TJCHYF1550	15	50	Complex				
	TJCHYF1630	16	30	Complex				
	TJCHYF1645	16	45	Complex				
	TJCHYF1840	18	40	Complex				
	TJCHYF1850	18	50	Complex				
	TJCHYF2040	20	40	Complex				

Beijing Taijieweiye Technology Co., Ltd


Add: 21 West Panlong St., Building 4, Mafang
Town, Pinggu District, Beijing, China 101204
Tel: 010-89471461 Fax: 010-89471460



Perfiller®

Expansible Embolic Coil System

CE 2292

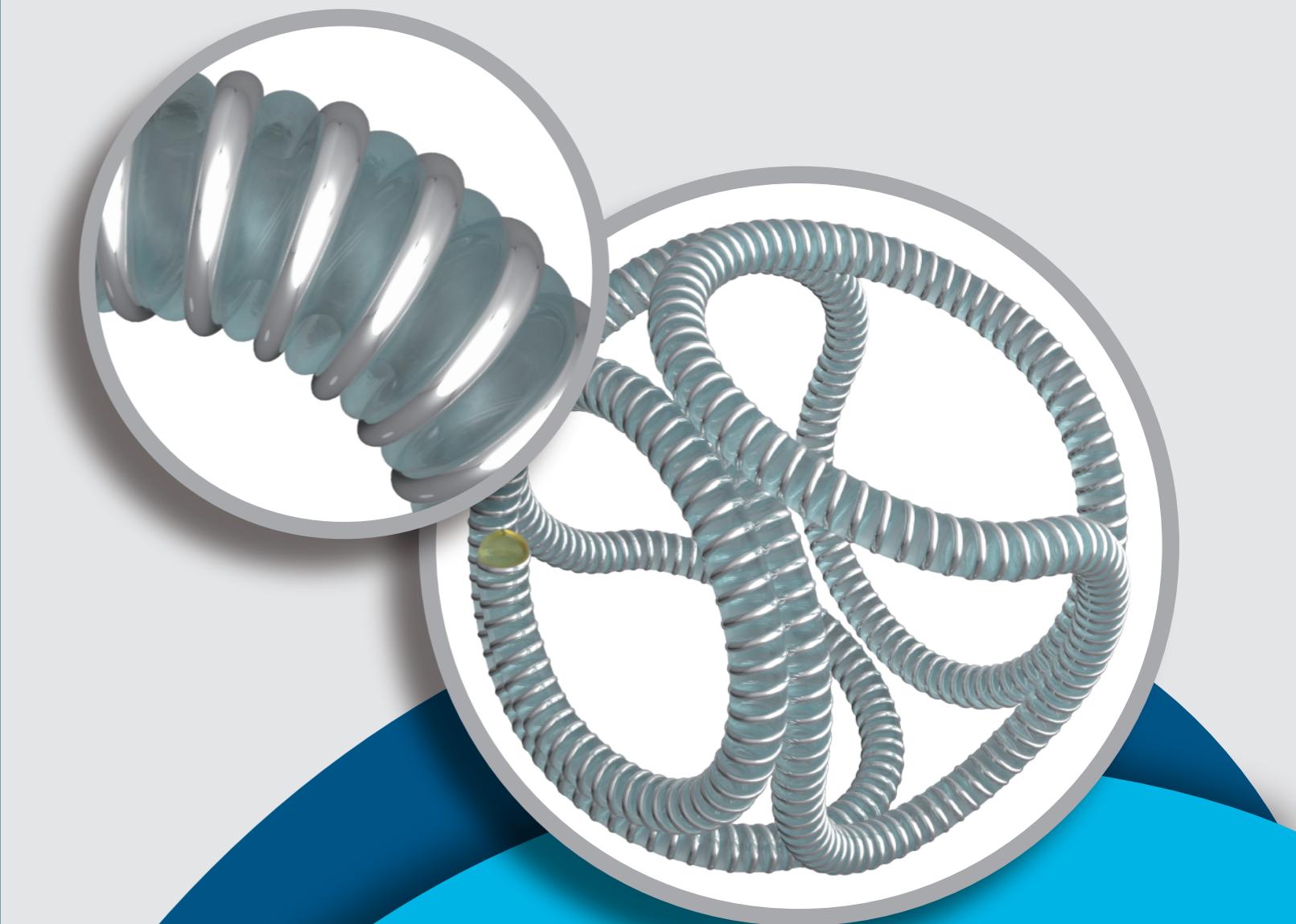

• No extra hydration process

• No limit on operation time