

This product is applied for the intracranial aneurysms, arteriovenous fistulas and other neurovascular abnormalities to be treated by interventional embolization.



Product Specification

Catalogue No. Helical (2D)	Catalogue No. Complex (3D)	Loop Dia. (mm)	Length (cm)	Catalogue No. Helical (2D)	Catalogue No. Complex (3D)	Loop Dia. (mm)	Length (cm)
TJCST1.502-2D	TJCST1.502-3D	1.5	2	TJCST4.506-2D	TJCST4.506-3D	4.5	6
TJCST1.503-2D	TJCST1.503-3D	1.5	3	TJCST4.508-2D	TJCST4.508-3D	4.5	8
TJCST1.504-2D	TJCST1.504-3D	1.5	4	TJCST4.510-2D	TJCST4.510-3D	4.5	10
TJCST0201-2D	TJCST0201-3D	2	1	TJCST4.512-2D	TJCST4.512-3D	4.5	12
TJCST0202-2D	TJCST0202-3D	2	2	TJCST4.515-2D	TJCST4.515-3D	4.5	15
TJCST0203-2D	TJCST0203-3D	2	3	TJCST0509-2D	TJCST0509-3D	5	9
TJCST0204-2D	TJCST0204-3D	2	4	TJCST0510-2D	TJCST0510-3D	5	10
TJCST0206-2D	TJCST0206-3D	2	6	TJCST0515-2D	TJCST0515-3D	5	15
TJCST0208-2D	TJCST0208-3D	2	8	TJCST0520-2D	TJCST0520-3D	5	20
TJCST2.502-2D	TJCST2.502-3D	2.5	2	TJCST0610-2D	TJCST0610-3D	6	10
TJCST2.504-2D	TJCST2.504-3D	2.5	4	TJCST0611-2D	TJCST0611-3D	6	11
TJCST2.506-2D	TJCST2.506-3D	2.5	6	TJCST0615-2D	TJCST0615-3D	6	15
TJCST2.508-2D	TJCST2.508-3D	2.5	8	TJCST0620-2D	TJCST0620-3D	6	20
TJCST0304-2D	TJCST0304-3D	3	4	TJCST0715-2D	TJCST0715-3D	7	15
TJCST0306-2D	TJCST0306-3D	3	6	TJCST0720-2D	TJCST0720-3D	7	20
TJCST0308-2D	TJCST0308-3D	3	8	TJCST0730-2D	TJCST0730-3D	7	30
TJCST0310-2D	TJCST0310-3D	3	10	TJCST0815-2D	TJCST0815-3D	8	15
TJCST0312-2D	TJCST0312-3D	3	12	TJCST0820-2D	TJCST0820-3D	8	20
TJCST3.506-2D	TJCST3.506-3D	3.5	6	TJCST0830-2D	TJCST0830-3D	8	30
TJCST3.508-2D	TJCST3.508-3D	3.5	8	TJCST0920-2D	TJCST0920-3D	9	20
TJCST3.510-2D	TJCST3.510-3D	3.5	10	TJCST0930-2D	TJCST0930-3D	9	30
TJCST3.512-2D	TJCST3.512-3D	3.5	12	TJCST1020-2D	TJCST1020-3D	10	20
TJCST0404-2D	TJCST0404-3D	4	4	TJCST1030-2D	TJCST1030-3D	10	30
TJCST0406-2D	TJCST0406-3D	4	6	TJCST1130-2D	TJCST1130-3D	11	30
TJCST0408-2D	TJCST0408-3D	4	8	TJCST1230-2D	TJCST1230-3D	12	30
TJCST0410-2D	TJCST0410-3D	4	10	TJCST1330-2D	TJCST1330-3D	13	30
				TJCST1430-2D	TJCST1430-3D	14	30
				TJCST1530-2D	TJCST1530-3D	15	30
				TJCST1630-2D	TJCST1630-3D	16	30
				TJCST1830-2D	TJCST1830-3D	18	30
				TJCST2030-2D	TJCST2030-3D	20	30

* Minimal microcatheter inner diameter compatibility: 0.0165"

EC	REP
----	-----

Lepu Medical (Europe) Cooperatief U.A.
Abe Lenstra Boulevard 36
8448 JB
Heerenveen
the Netherlands


Beijing Taijiweiyi Technology Co., Ltd

Add: 21 West Panlong St., Building 4, Mafang Town,
Pinggu District, Beijing, China 101204
Tel: 010-89471461
Fax: 010-89471460

Please read the Instruction for Use for indications,
contraindications and possible complications.

Perdenser®

Embolic Coil System

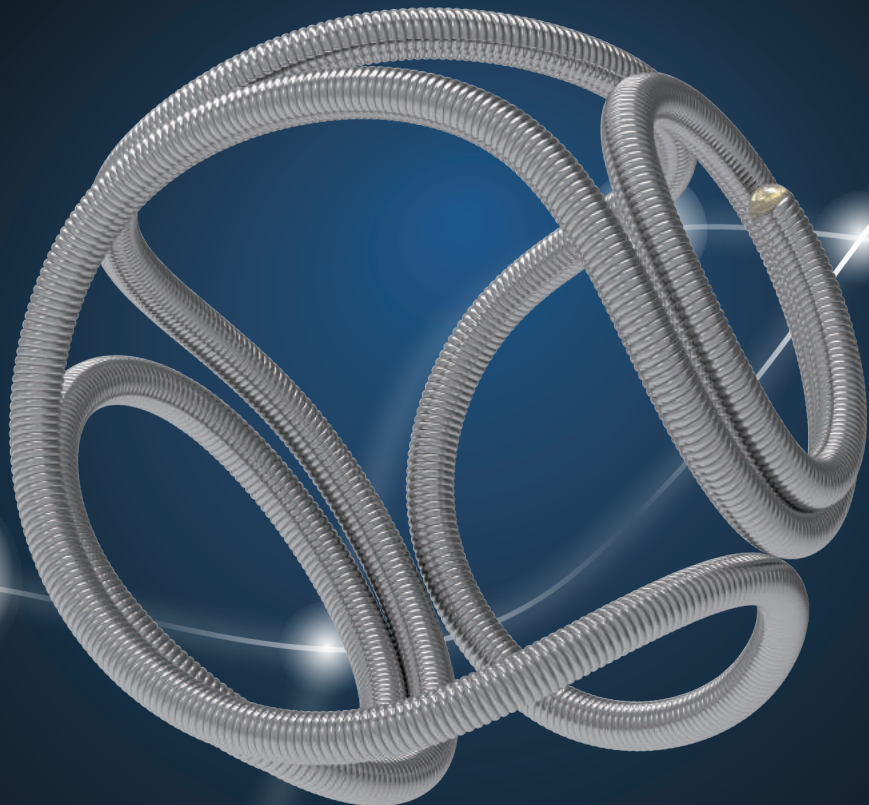

LEADING TECHNOLOGY • CHERISHING LIFE

1/6 Full Range of Sizes for Accurate Positioning

Full range of sizes of coils with both helical and complex structures. The coil system is compatible with 0.0165"-0.021" microcatheter inner diameter. Two markers design provides accurate positioning for both delivery system and detaching point.



Framing



Filling




Closing




3/6 Strengthen Anti-unwinding Force

The maximum anti-unwinding force could reach up to 0.5N that provided by the polymer anti-unwinding system, which guarantees its performance for coil readjustment.

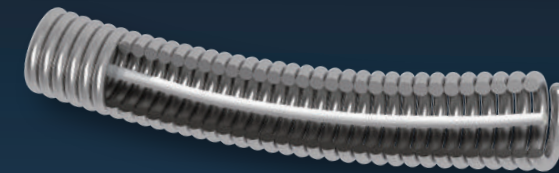


4/6 Soft Shrunken Distal Loops

Pt-W alloy coil and round polymer tip enhance its softness. The first loop of the coil is shrunken to minimize the risk of coil prolapsing.


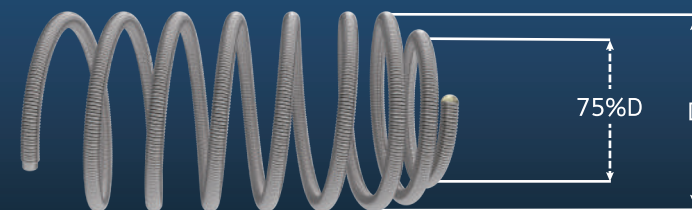


2/6 Multi-segments of Delivery System

The delivery system is designed with gradient hardness, which provides stiff proximal part and soft distal part. The longer soft segment combines with shorter detaching segment enhances the flexibility of the delivery system to go through the tortuous vessels smoothly. The stiff proximal shaft provides better pushability for easier manipulation.

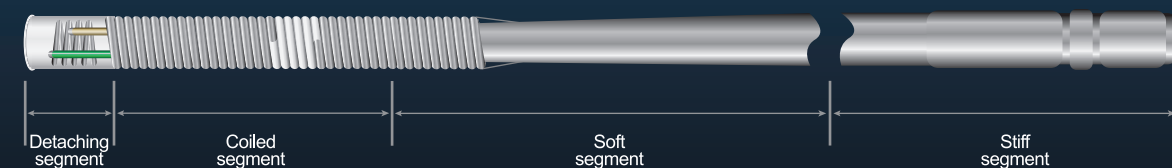

5/6 Uniform and Compact Filling Performance

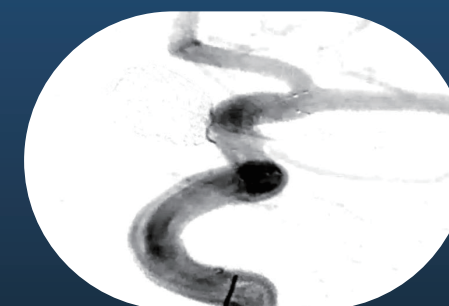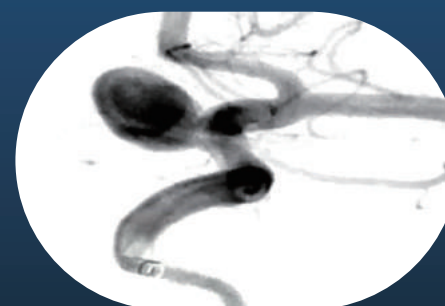

Optimized Ω open loop design provides random yielding points to fit any shape of aneurysm. The coil spreads out to any direction to form a stable frame, and advances centripetally to full fill the space.

6/6 Safe and Reliable Immediate Detaching System

The detaching point is applied with polymer wire and immediate fusing technique, which ensure the detachment finishing within 1 second. The distal end of the coil is smooth and soft after detachment without any additional chemical residues. The operation of the detachment is easy, stable and safe.

